


# Glass Water distiller

## LGWD-A10

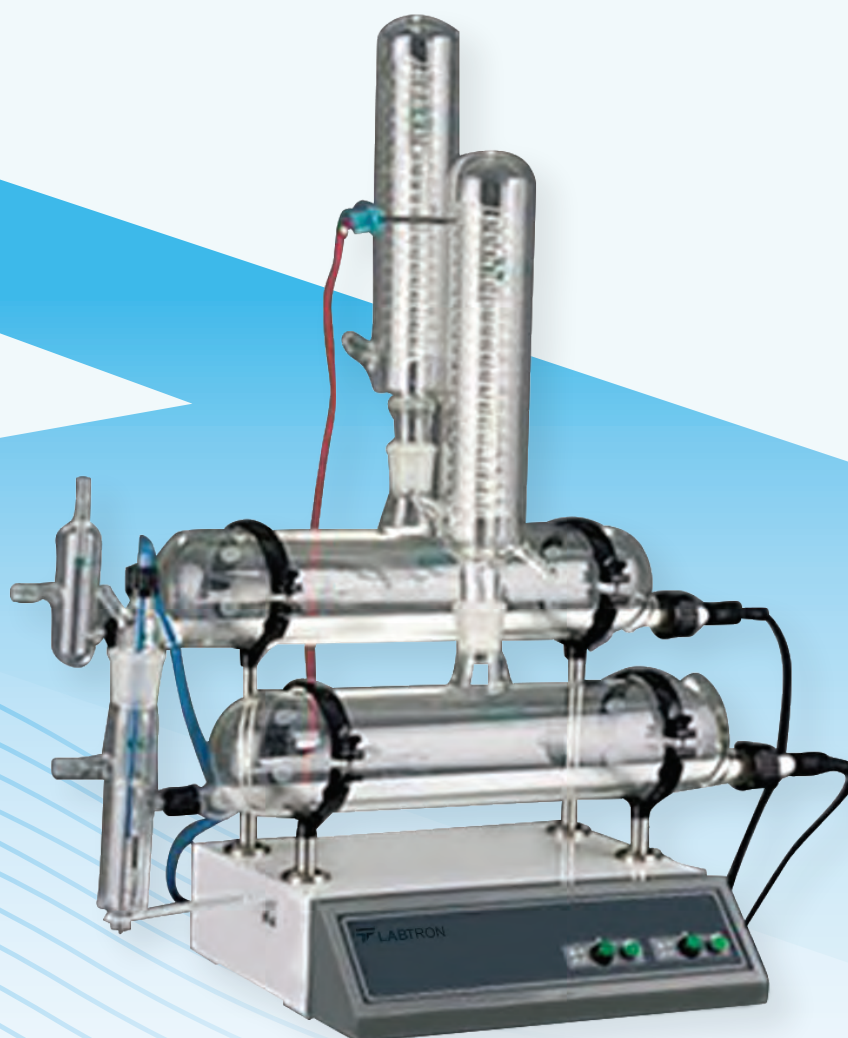

## Glass Water distiller LGWD-A10

Glass Water distiller LGWD-A10 is a dual water distiller, equipped with quartz glass heater for saving energy, for easy cleaning and long life. High purity water is obtained through the secondary distillation. High quality glass distillers are used to avoid contact between water and metal part. Distillation switch is controlled through temperature controller and a reed water leveler.

### Features

- Dual water distiller
- Quartz glass tube heating type
- Relay protection device for safe operations
- High quality wear resistant SS Stainless steel support column, column base
- High quality glass condensation cycle system
- Automatic switch on/off mode to maintain water level height at safe operations

### Application

Used to obtain pure distilled water required for different preparations across research and testing laboratories of pharmacy, biotechnology, tissue culture, microbiology, chemistry and other fields.

### Specifications

Model no.	LGWD-A10
Water discharge	1600 ml/ h
Input power	3.0 kW
Outer dimension	700 × 340 × 800 mm
Power supply	AC 220 V, 50 Hz
Net weight	12 kg
Gross weight	14 kg