


Paraffin Trimmer LPFT-A1 Series

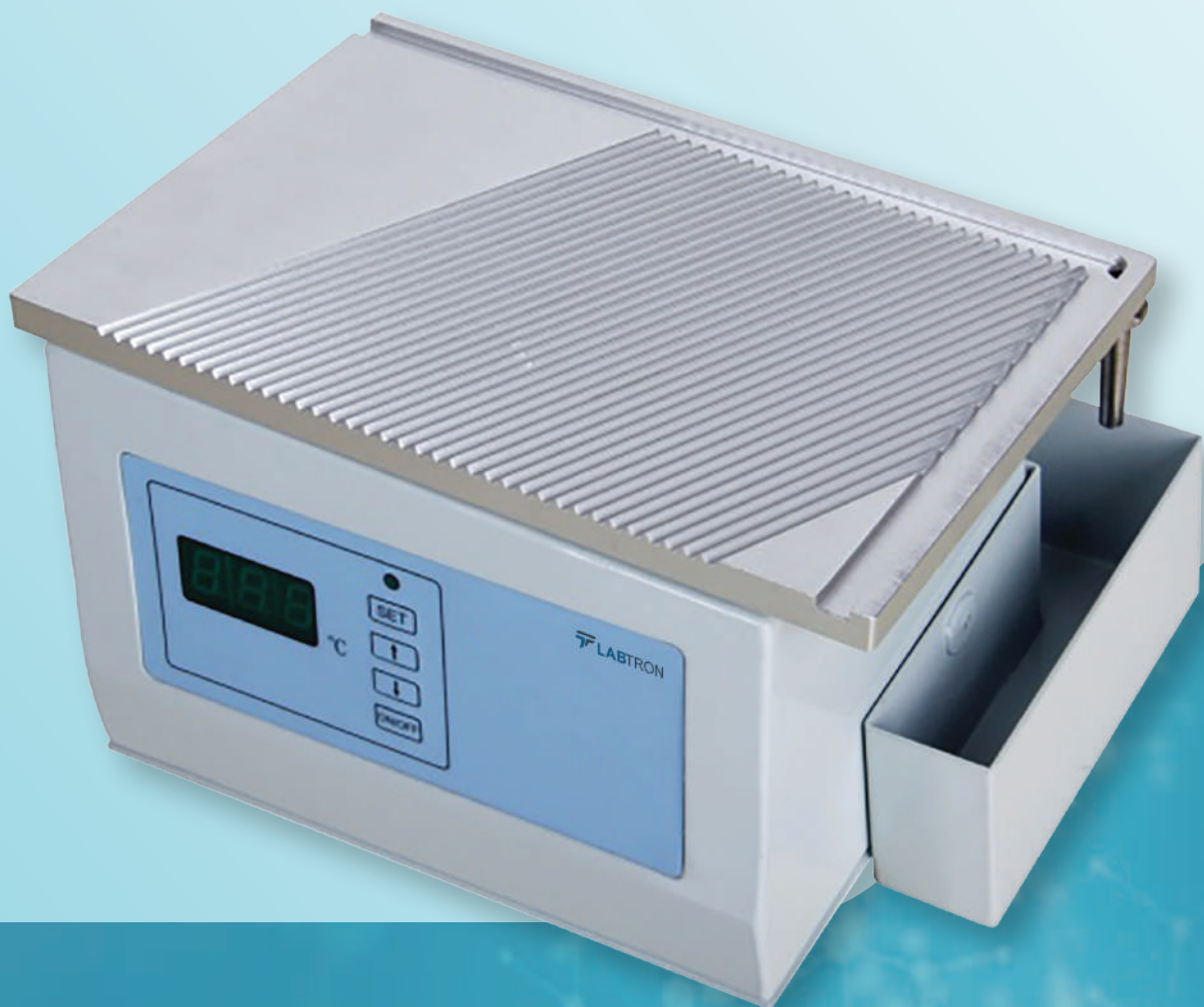

Paraffin Trimmer LPFT-A1 Series

Paraffin Trimmer LPFT-A10

Paraffin Trimmer LPFT-A10 is a compact tabletop trimmer for removing excess paraffin from embedding tissue cassettes. Melted paraffin is collected into the attached waste collection box. Trims multiple blocks simultaneously. Equipped with grooved surface for proper wax drainage. It is ideal used for removing excess paraffin from cassettes in the medical, histology, research laboratories.

Features: (LPFT-A10)

- Microcomputer controller
- Temperature sensors are made precise and reliable
- Automatically save the last operational settings
- Memory and automatic restoration function
- LED display screen for temperature
- Equipped with melted wax collection box
- Grooved surface for proper wax drainage

Applications

Widely used in removing excess paraffin from tissue, embedding, biopsy, mega cassettes in the medical institutions, research labs, hospitals, pathology laboratories, etc.

Specifications

Model no	LPFT-A10
Temperature range	RT to 90 °C
Temperature accuracy	± 1 °C
Working area	230 mm × 160 mm
Power consumption	150 W
Voltage	AC 110 V/ 220V ± 10 %, 50 /60 Hz
External dimension (W × D × H)	255 × 180 × 155 mm
Packing dimension (W × D × H)	330 × 260 × 230 mm
Weight	3 kg

Standard Accessories

Model no	Accessory name	Quantity
1.	Waste paraffin collection boxes	5